Problem List

Last Modified on 10/22/2024 12:12 pm EDT

The Problem List chart tab is used to capture a list of conditions a patient has been diagnosed with. This list of problems is used to facilitate workflows throughout the application including the Superbill and Order Entry chart tabs. The Problem List chart tab will also double as the assessment for many users by adding notes to a problem via the edit window.

Problem List Overview

Problem Search

The **Add Problem** search bar allows a user to find a diagnosis that you would like to add to the selected patient's chart. When a user begins to type in the search box, a list of results are displayed.



The list of returned results are categorized as either ICD, Problem or Full.

• ICD: displays a list of codes that have been added to your practice via the Diagnosis Code Setup window.

Please note, the list of Diagnosis Codes in your database are customized for your practice. You may find it necessary to add, change, or remove diagnosis codes to meet the unique needs of your practice.

- **Problem**: displays a list of SNOMED CT codes that are specifically designated as common problems. This is a helpful search tool because SNOMED provides "physician friendly" names and also supports searching by abbreviations (e.g. CKD). SNOMED CT even cross-walks to ICD codes.
- Full: displays a full list of SNOMED CT codes. These codes may or may not cross-walk to an ICD code.

User Settings

Search List Type

Favorite Problem Search

If a user selects the drop down to the right of the search box, a list of favorite ICD codes will be displayed. This list of favorite codes is generated by looking at all problems that were added across the entire patient population for a given time frame. There are settings that you can change to control this experience, see the User Settings list below.

Company Settings

Auto codify Problem List Codes to SNOMED

User Settings

- Problem List Favorite Search Type
- Problem List Favorite Search Days
- Problem List Favorite Search Max Rows

Problem List Summary

The **Problem List Summary** contains diagnoses that have been added to the patient's chart and are currently in an **Active** status. You can interact with these problems by clicking on them.

D291 Benign neoplasm of prostate (disorder)

06/18/2017

New sx occured today. Notice odd pain in left leg that was not there before. Will send for MRI and recheck in 3 weeks.

08/20/2019 5:51PM

The Problem List Summary will display the following elements for each code:


- ICD Code: This is the ICD-9 or ICD-10 code associated with the problem. Not all problems have a corresponding ICD code.
- Problem Description: This is the name of the selected condition.
- Onset Date: This is the date that the problem began or was diagnosed.
- **Notes:** On the summary view, you will see the last note added for this problem. There can be one to many notes on any given problem.
- Last Modified Date: This is the last time this problem was changed.
- Priority Level: A number representing the priority order for the condition is displayed.
- Interaction Icon: You may see a yellow or red triangle next to a condition. This indicates an interaction between this problem and a patients medication list.

User Settings

Problem List Summary View

Problem List Edit Window

The Problem List Edit window is accessed when a new problem is added or a problem is selected from the Problem List Summary. This allows a user to access and update the details related to an individual Problem List item.

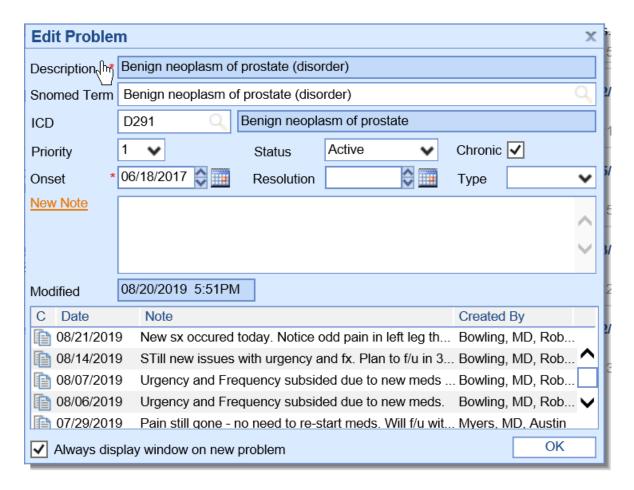

The edit window will display the following elements:

- **Description**: This is the description of the problem as seen on the screen. This is pulled from the database the code was selected from.
- **Snomed Term:** A problem may or may not be linked to a SNOMED term. This typically happens automatically when selecting a problem from the 'Problem' or 'Full' option within the search.
- ICD: This is an ICD code lookup. This is searching the Diagnosis Code Setup table.
- ICD Description: This is the corresponding description for the selected ICD code.
- Priority: This is the priority for the given problem.
- Status: This is the current status for the problem.
 - o Active: The problem is actively being treated and will display on the summary.
 - o Inactive: The problem is not active and will not display on the summary.
 - Removed: The problem is no longer being treated and will not display on the summary.
 - Resolved: The problem is no longer being treated because it has been resolved. This will not display on the summary.
- Chronic Flag: This indicates that the condition is chronic or acute.
- **Onset**: This is the date the problem began. This will default to the current date when a new problem is added.
- **Resolution**: This is the date the problem was resolved. Setting this date will auto-set the status to 'Resolved'.
- **Type:** This is used to define a problem as the patient's cognitive or functional status. This is typically left blank.
 - o Cognitive
 - Functional
- New Note: This allows a user to type or dictate a new note for the problem. This note will be time and
 date stamped. If a new note is added on a problem and that problem is used on the superbill that day, this

note will be added into the SOAP note.

- Modified: This is the date/time that the problem was last changed.
- **Note List:** This is a list of all notes that have been added to the problem over time. Clicking on a note will allow you to see a more detailed view. You can even copy forward a previous note into today's note.
- Always Display Window on New Problem: When selected, this will ensure that the Edit Window display each time a new problem is added from the search. If deselected, the edit window will not appear.

Problem List History


To see a complete list of all problems for a patient, regardless of status, a user can click the **History** button on the toolbar.

History						
Alert	Status	Code	Description	Onset	Modified	Note
	Active	D291	Benign neoplasm of prostate (06/18/2017	08/20/2019 5:51PM	New sx occured today
	Active	N401	Benign prostatic hyperplasia (08/22/2018	08/22/2018 11:18AM	,
\triangle	Active	E1121	Diabetic renal disease (disord	11/05/2018	11/05/2018 1:53PM	Being treated by PCP D
	Active	N200	Kidney stone (disorder)	11/13/2018	07/29/2019 2:21PM	SCheduled an ESWL
	Active	N3281	Bladder muscle dysfunction	02/12/2019	02/12/2019 1:37PM	Worse after subdural h
\triangle	Resolved	R310	Gross hematuria	11/06/2018	02/12/2019 1:13PM	Resolved
\triangle	Resolved	R31	Hematuria	06/18/2017	02/12/2019 1:13PM	Resolved
	Inactive	N20	Calculus of kidney and ureter	06/18/2017	10/02/2018 12:15PM	Inactivated
	Inactive	R3919	Slowing of urinary stream (fin	06/18/2017	04/03/2018 5:08PM	Inactivated
	Inactive	R972	Raised prostate specific antig	06/18/2017	07/06/2018 2:30PM	Inactivated
	Inactive	N400	Large prostate (finding)	07/03/2017	02/08/2019 9:35AM	Inactivated
	Inactive	N451	Acute epididymitis (disorder)	10/12/2017	01/16/2018 10:54AM	Inactivated
	Removed	R9720	Elevated prostate specific anti	05/13/2018	06/06/2018 5:50PM	Removed


Nearly all of the same features and functionality that apply to the **Problem List Summary** window apply to the history. The major difference is that the history will display all problems including those that have been removed, resolved, and inactivated.

Add a New Problem


1. Once in a patient's chart, select the **Problem List** chart tab.

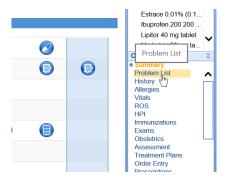

2. Begin typing in the Add Problem search box or select the drop down to display the Favorites List.

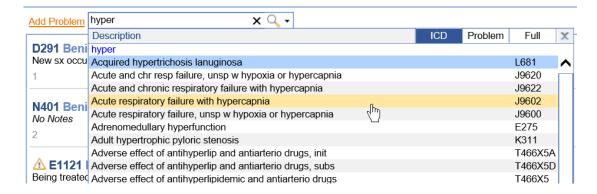

3. Once an item is selected from the list, the Edit Problem window will display. Add or change details as needed.

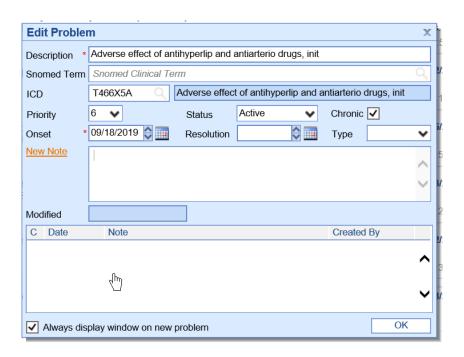

4. Click **OK** when you are done changing this problem.




5. Press **Save** when you are done adding all other problems.



Remove a Problem

To remove an active problem from a patient's chart, the status of the problem must be changed. Problems that are visible in the **Summary** view of the Problem List have a status of Active. If a problem has any other status, it is displayed in the **History** view of that chart tab. Follow these steps to remove a problem from the active list:


1. Once in a patient's chart, select the **Problem List** chart tab.

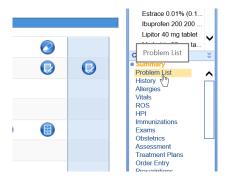

2. Click on the problem item from the grid on the **Summary** window. An **edit window** will appear.



3. Change the Status drop down to anything other than Active.

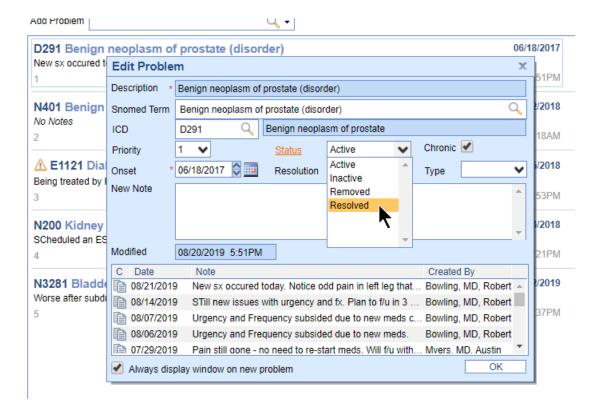

4. Click the Save button




Create an Assessment using the Problem List

The Problem List has the ability to be used to create an Assessment/Plan for a specific condition. To do this, a user will make a comment on a problem. The user will then choose to use this diagnosis on their Superbill for that same day. When the comment and the Superbill are created on the same day, the comment will display on the note.

Step 1: Add the Assessment to the Problem


1. Once in a patient's chart, select the **Problem List** chart tab.

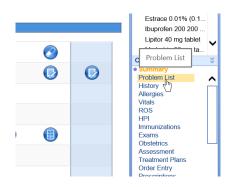


2. Click on the problem item from the grid on the Summary window. An Edit Problem window will appear.



3. In the New Note field, type, dictate, or copy forward an old note.

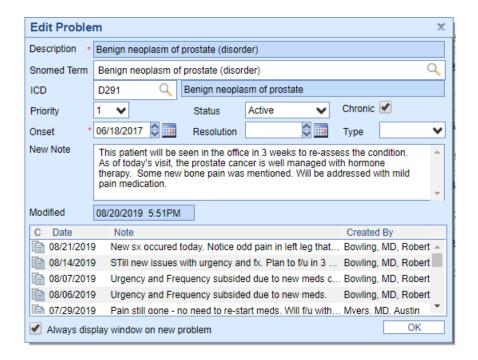

4. Click **OK** on the Edit Problem window.

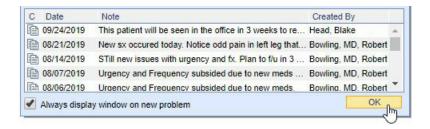


5. Click Save on the Problem List chart tab.




Step 2: Add the Problem to the Superbill


1. Select the Superbill chart tab

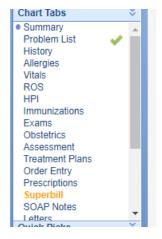


- 2. Add the problem to the **Charge Summary**. This can be done in many different ways.
 - 1. If you created the assessment on the problem today, the system should recommend that you add the Diagnosis to the Superbill for this patient.

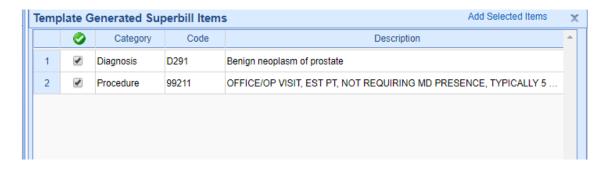

You can manually select the diagnosis from the **Problem List** bluebar that is automatically generated.

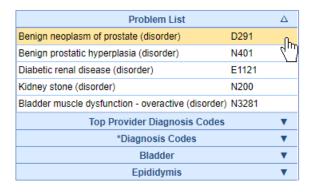

3. You can search for the diagnosis in the **Lookup** bluebar.

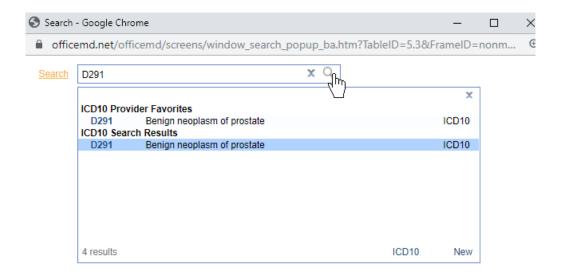


3. Add other necessary details to the Superbill. Ensure the problem that you performed the Assessment on is listed in the **Charge Summary**. Save the Superbill.



Step 3: Create the Note

1. Select the **SOAP Note** chart tab

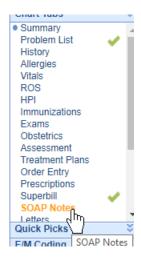

2. The note should auto-populate. If not, click Re-Populate Note.



3. You will notice the note that was added to the problem will show in the **Assessment** portion of the note by default.

Assessment

Diagnosis

D291: Benign neoplasm of prostate (This patient will be seen in the office in 3 weeks to re-assess the
condition. As of today's visit, the prostate cancer is well managed with hormone therapy. Some new
bone pain was mentioned. Will be addressed with mild pain medication.)

Procedures


 99242:OFFICE CONSULTATION, 3 KEY COMPONENTS:EXPAND PROB FOCUS HX;EXPAND PROB FOCUS EXAM;STRTFWD MED DECISN

Plan

View the Problem List History

The **Problem List History** screen displays a patient's complete list of diagnoses regardless of their status (Active, Inactive, Removed or Resolved). Follow these steps to view the Problem List History:


1. Once in a patient's chart, select the **Problem List** chart tab.

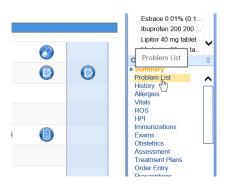


2. Click on the History button on the EMR Toolbar.




3. Click on a problem to open the Edit Problem screen.

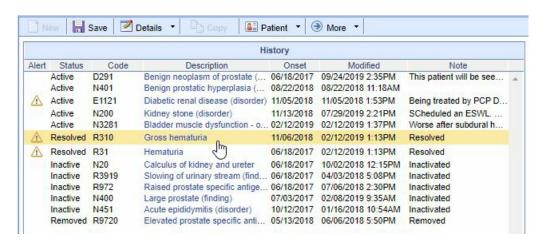

4. Edit the problem as necessary.

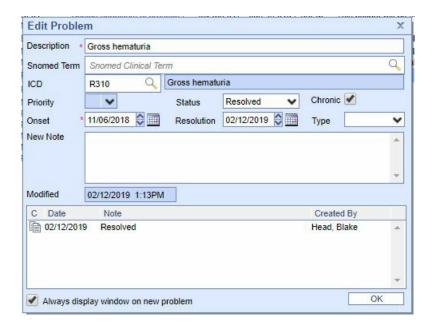


5. Press the **OK** button.

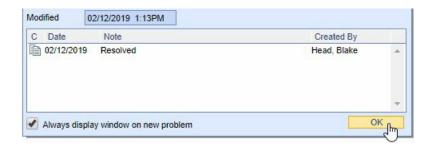

6. Click **Save** when you are done editing/viewing all other problems.



7. Click the **Details** button to return to the **Problem List Summary** screen

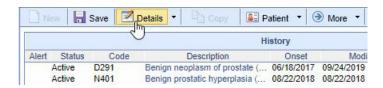


User Problem List Favorites

COMING SOON (Release: November 1, 2024)

Introducing the new Favorites Feature in the Problem List!

Users will now notice a star icon next to each problem, allowing them to mark items as personal favorites. This star makes it easy to track and quickly access frequently used problems.

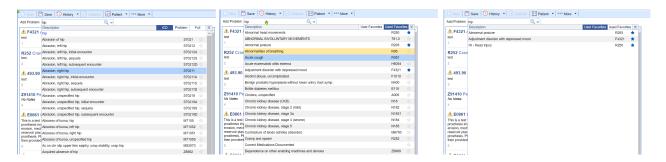

From the search bar or by using the drop-down menu in the search bar, users can access their saved lists and select or deselect the star icon to manage favorites. Items with a grayed-out star can be marked as User Favorites* by clicking the icon, turning it **BLUE**. These items will then appear in the User Favorites tab.

If a user is in the User Favorites tab and clicks the star icon next to an item, the item will be deselected, the star will turn to a gray outline only, and it will no longer appear in the User Favorites tab.

781.0

R06

H6054

F4321

F1010


E119

A009

N182

N184

N1831

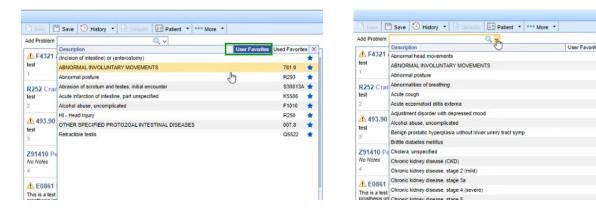

By default, the Used Favorites tab is displayed, but users can set the User Favorites tab as their default view through a new user settings option (*User Setting:* Default favorite list to User Favorites). This makes managing and accessing personalized lists even easier!

